


Labycare-Tube Stripper. (SKSS388B)

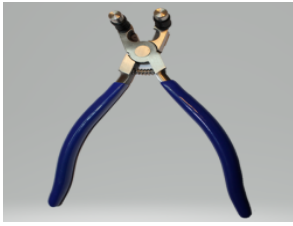

Donor Tube Stripper/Sealer is used to strip blood from donor tubing into the blood bag, seals the tube with clip and cuts tubing in the preparation of blood segments from the blood collection and processing procedures. The blood in the tubing segments is stripped into collection bags. Using the blades, the segments are further separated by cutting between the seals. The body is made of stainless steel with Teflon rollers.

Description :

Tube Stripper, easy to use and ergonomic, can be used with different sizes of Blood Bag tubes due to its screw that allows adjusting the distance between the rollers.

The body is made of aluminum with stainless steel rollers. Its lightweight body makes it extremely comfortable for everyday use and resistant to falls.

The stripper includes the spring for automatic return slot and a blade to cut the tubes.

Specifications :

- High-quality stainless steel
- Anti-rust, lightweight and sturdy
- Ease of cleaning by wiping or soaking in IPA
- Ergonomic design for good grip and stripping
- Ease fatigue over prolong use
- Self-centering Teflon roller system keeps the tubing in position
- Smooth traction for effective stripping of blood components
- Internal spring system keeps handles in open position
- Minimize entanglement with blood tubing
- Include 50 Terumu aluminum seals

